

Valuable and Accessible IT Solutions


Mindray ePM devices can connect to the Central Monitoring Station (CMS) and eGateway through both wired and wireless networks, as well as interfacing with third-party electronic medical records (EMR) via HL7 output directly.

The ePM helps enhance clinical work flow and efficiency with it's flexible yet reliable connectivity capabilities.

- The View Other Patient function allows caregivers to see, in real-time, up to 12 other beds on a single ePM screen. This seamless information exchange between bedside monitors can help caregivers view all their patients at once, without the need for CMS.
- With the ePM Caregroups function caregivers can quickly find and review their assigned patients or ward when connected to the CMS.

Data from ePM devices can be easily connected to the CMS and Mobile Viewer, giving clinicians access to their patient data anytime and anywhere in the hospital.

- The CMS Early Warning Scores (EWS) dashboard provides an intuitive display of patient status, with dynamic updates pushed to the Mobile Viewer, alerting caregivers to changes in patient conditions and potential risk of deterioration.

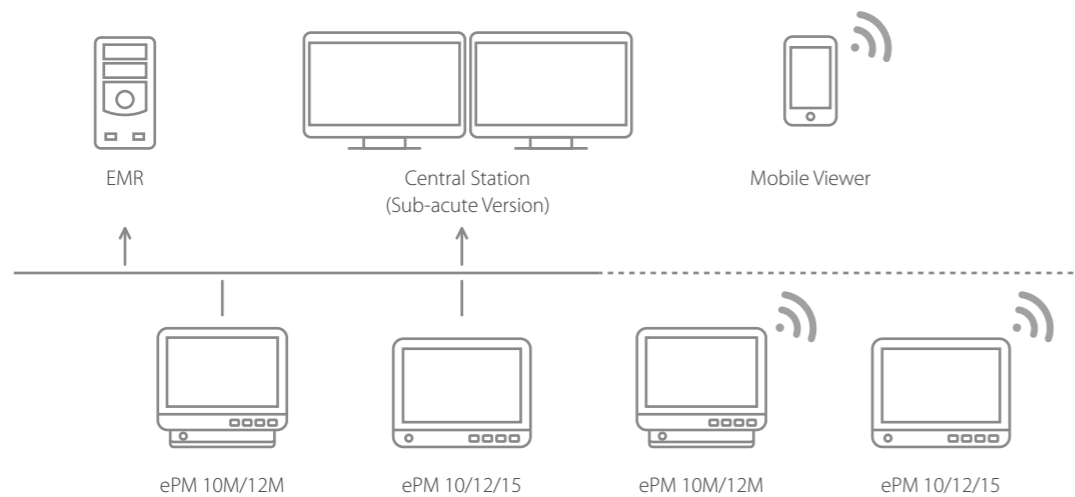

ePM 10M/12M

Modular Patient Monitor

The evolution of simplicity



MRB015V221




Inspired by the needs of customers, Mindray patient monitors adopt advanced technologies and transform them into accessible innovation. The ePM delivers excellent visual experience, intelligent operation, accurate physiological measurements, smooth workflow and comprehensive connectivity options for demanding hospital settings, such as Emergency Rooms, Recovery Units, Sub-acute Units and General Wards.



Minimalist Design



Multi-touch capacitive screen
Supports gestures operations

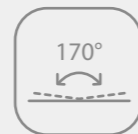

Wide viewing angle
Makes display more visible



1280x800 pixel (10.1"/12.1")
Provides HD visual experience



Auto brightness
Reduces light interference at night



Fanless design
Reduces the risk of cross-contamination



Durable and robust casing
Validated for cleaning with 49 leading disinfectants

Extensive Data Storage at the Bedside

2400 hrs
Trends @ 10 minutes

48 hrs
Full disclosure waveforms for all parameters

3000 sets
NIBP measurement

2000 sets
Events

400 sets
OxyCRG events

Note: These are the Maximum storage capacity of ePM devices with 16G storage.

Expand Parameters to Meet Your Clinical Needs

- Expanding extension rack supports plug-and-play parameter modules, which are also compatible with Mindray BeneVision N series monitors
- Ergonomic, concealed handle without cleaning blind spot




Wireless Monitoring

- ePM 10/12M can upgrade the receiver module to wearable sensor for wireless monitoring, so as to solve the inconveniences caused by monitoring cables and make patients more comfortable.
- Unnoticeable wearing experience
- Reliable design. 1.5m Drop protection. Resistant to 48 kinds of disinfectants. 36h Long battery life
- Medical-grade multi-parameter monitoring
- Innovate health parameter monitoring. Monitoring patient's exercise time, sleep time and fall-down detection



Medical-grade multi-parameter monitoring

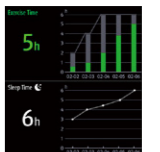

Health parameter monitoring

Flexible Mounting Solutions

- A wide range of mounting solutions designed for various clinical settings
- The release mechanism allows for quick removal from the wall mount or rolling stand for transport



Accurate, Reliable Parameters


Comprehensive Monitoring

Integrated Platinum™ MPM parameters:
3/5/12-lead ECG, respiration, SpO₂*, temperature and NIBP.

- Multi-lead(4) ECG algorithm with ST & QT analysis
- Low perfusion SpO₂ algorithm
- Fast, accurate and motion tolerant NIBP algorithm, validated by British Hypertension Society (BHS)

Wide measurement range and anti-interference performance ensures excellent parameter accuracy and reliability.

ePM also provides advanced parameter options:
4-ch invasive blood pressure, EtCO₂ and cardiac output, multi-gas and BIS, making it suitable for a wide range of clinical settings.

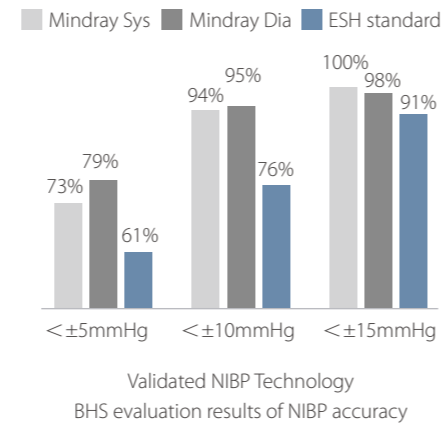


* Mindray provides 3 options of SpO₂ measurement, Mindray, Masimo and Nellcor. For further information about the availability of Masimo and Nellcor SpO₂, please contact with your local sales representatives.

Reduce False Alarms with CrozFusion™



Innovative multi-parameter alarm analysis can reduce false arrhythmia alarms and promote the accuracy of heart rate and pulse rate, and help to alleviate alarm fatigue.

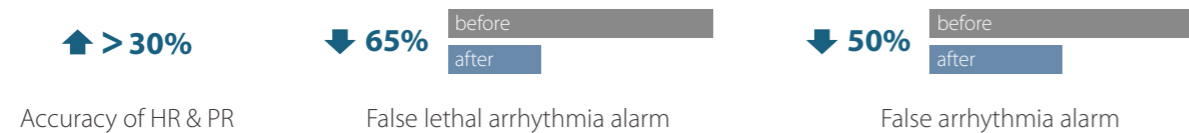

Note: The results are based on an evaluation by Mindray multi-parameter fusion database.

Intelligent Alarm Strategy

- Alarm highlight: Using special fatal alarm sound. Highlight the fatal alarm
- Alarm limits recommendation: Support personalised threshold settings for different patients based on trend data
- ARR Alarm refractory period: Avoid repeated meaningless alarm



Simplicity at Your Fingertips

Intelligent Operation Experience

- Operate with gestures, just like a tablet PC
- Access to the most common functions in 2 steps or fewer
- Quickly identify disconnected sensors with the innovative AlarmSight technology
- Screen lock for easy cleaning



Smooth Workflow

Based on clinical insight, the ePM has optimised workflows to support caregivers at the bedside, swiping the touchscreen to switch between commonly used functions and interfaces, enabling clinical tasks to be completed quickly and accurately.



View from a distance
Intuitive big numerics



View at bedside
Highlights abnormal readings



Ward rounds or nurses hand over
Quickly review the patient status changes




Review and analyse
24hrs waveform review and critical alarms

Early Warning Scoring (EWS)

Mindray ePM monitors provide a point-of-care EWS calculator to help clinicians track and document signs of patient deterioration, aiding faster and more informed patient care decisions.


- Compliant with the National Early Warning Score (NEWS), National Early Warning Score 2 (NEWS2) and Modified Early Warning Score (MEWS) protocols
- Alternatively, create custom protocols to suit your hospital requirements
- Intuitive layout and trends review helps caregivers visualise data quicker
- Automate EWS calculations quickly at the bedside
- Display score escalation instructions on-screen to remind caregivers make rapid care decisions
- Integration to the Electronic Medical Record (EMR) for fast, accurate electronic vitals and early warning scoring documentation

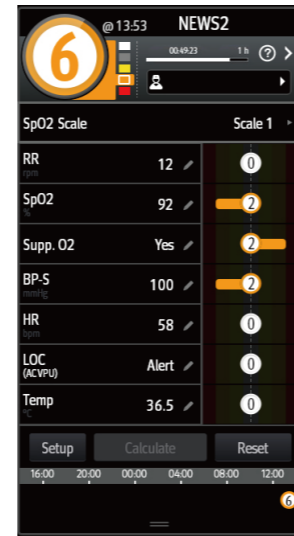

The Early Warning Score tool, as displayed on ePM devices

Clinical Assistive Applications (CAA)

The ePM provides efficient Clinical Assistive Applications (CAA) to help support safe and efficient decision making in mid-acuity and general ward areas.

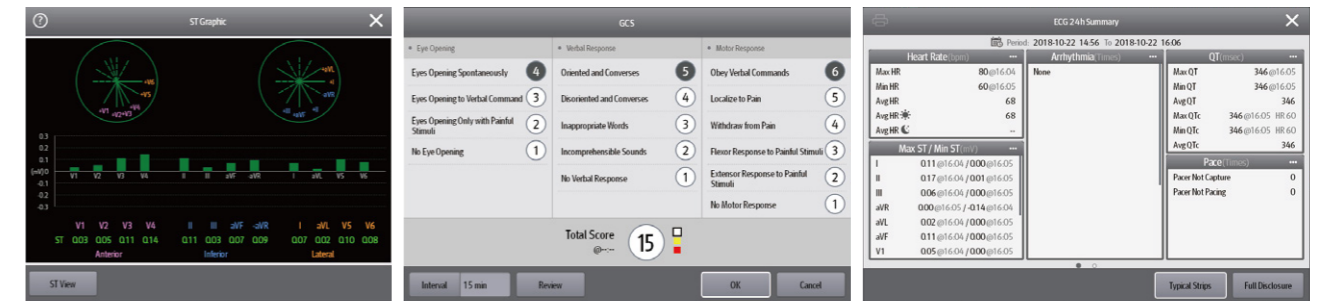

ST Graphic™

Glasgow Coma Scale

24 hours ECG summary

Supporting Safety in Neonates

Targeted Goal Screen

- Target parameter dashboard, visualised goal management such as oxygen therapy.
- 24 hours of target parameter statistics helps caregivers to evaluate the treatment effects
- Target parameter trend helps clinicians to identify patient situation. Highlighted part is combined with certain bar in the histogram



OxyCRG

- Effectively identify apnea of prematurity as ABD event
- Detailed and complete records of events help caregivers quickly identify the cause

